

DISCOVER POWERFUL WASHING

WITH VERSATILITY AND INCREASED EFFICIENCIES




**VISION®
1300**
Series Cart and Utensil
Washer/Disinfectors

Flexibility comes standard – from 3 unique models to an optional instrument cycle, STERIS Cart and Utensil Washers/Disinfectors help ensure your productivity goals are consistently met throughout the ebbs and flows of the day. Using significantly less water*, they minimize their impact on the environment, while making operation effortless with an intuitive design.



PRODUCTIVITY

- Provides additional capacity for reprocessing instrument trays and loaner sets. The Vision 1300 Series Cart and Utensil Washer's/Disinfector's optional instrument cycle features longer spray arms with nozzles at a 45 degree angle to ensure complete coverage

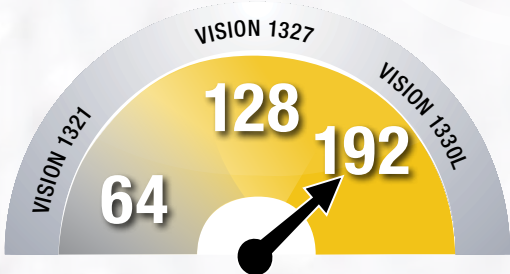

Carts per 8 hour day**



- Eliminates wait times. The Vision 1300 Series Cart and Utensil Washer's/Disinfector's HEPA-filtered drying system and automatic floor tilt deliver thorough drying immediately upon cycle completion



PERFORMANCE

- Our Full Coverage Spray Pattern is the ultimate cleaning technology with its cross pattern spray arm nozzles, impingement options, and dual speed pump. All of these elements work in concert to create our Full Coverage Spray Pattern, giving you the assurance you need for consistent performance.
- Provides thermal disinfection with each Standard Cycle



GREEN SAVINGS

Our Cart washers are built to run efficiently and minimize their impact on the environment. The Vision 1300 Series Cart and Utensil Washers/Disinfectors have been designed with water and energy-efficient features, providing the lowest water consumption per cycle of any washer on the market

- Full tempered glass door is made with 20% recycled glass
- 70% of water is recycled when cart cycle is used


LOWER TOTAL
COST OF
OWNERSHIP

Vision 1300 Cart and Utensil Washer/Disinfector uses up to

75% LESS WATER

per cycle*

CONCENTRATED AND
BIODEGRADABLE
CHEMISTRIES

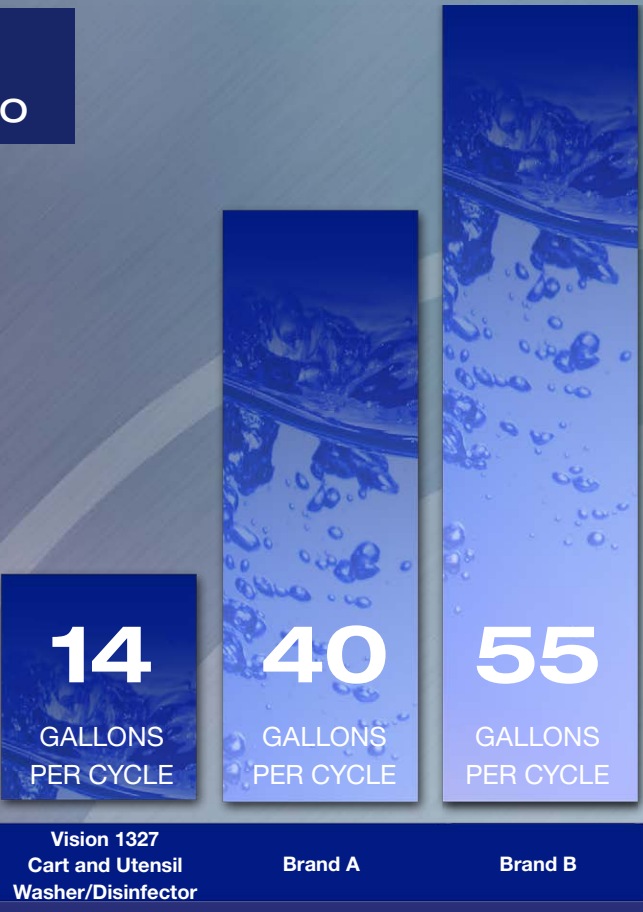



EASE OF USE

- Integrated cycle printer on the clean side
- Giving full visibility of carts or instruments, the Vision 1300 Series Cart and Utensil Washer's/Disinfector's automatic sliding glass door is easy to operate



- Increase staff confidence and safety with our intuitive LCD large touch-screen control panel which provides clear visibility and access to critical information



- Cool touch handles allow for safe and immediate handling of loads following completion of cycle
- Maintain compliance with our USB download port for easy access to cycle data

- The instrument rack is conveniently designed for easy loading and unloading from all sides
- Integrated flow meters provide verification of chemistry injection rates in our Vision 1300 Series Cart and Utensil Washers/Disinfectors
- Safety features include chamber illumination, emergency stop and door interlock

Size (W x H x D):

Vision 1321 Cart and Utensil Washer/Disinfector

Maximum Overall Dimensions:

120 x 112.75 x 106" (3048 x 2864 x 2692 mm)

Effective Chamber Load Capacity:

40 x 80 x 82" (1016 x 2032 x 2083 mm)

Vision 1327 Cart and Utensil Washer/Disinfector

Maximum Overall Dimensions:

120 x 112.75 x 130" (3048 x 2864 x 3302 mm)

Effective Chamber Load Capacity:

40 x 80 x 106" (1016 x 2032 x 2692 mm)

Vision 1330L Cart and Utensil Washer/Disinfector

Maximum Overall Dimensions:

120 x 112.75 x 142" (3048 x 2864 x 3607 mm)

Effective Chamber Load Capacity:

40 x 80 x 118.1" (1016 x 2032 x 3000 mm)

Vision 1300 Series Cart and Utensil Washers/Disinfectors are a part of an unbroken chain of innovative solutions from STERIS.

*Compared to non-STERIS Cart Washers/Disinfectors and previous STERIS Cart Washer/Disinfector models

** Number of carts based on Quick Cycles and size of cart

That is the power of...



For the washer/disinfector that fits your department's needs, contact your STERIS Representative or visit:

<http://www.steris.com/products/washing-and-decontamination-systems>

**ISO AAMI
15883
COMPLIANT**

STERIS

STERIS Corporation
5960 Heisley Road
Mentor, OH 44060-1834 · USA
440-354-2600 · 800-548-4873
www.steris.com

©2016 STERIS Corporation.
All rights reserved. Printed in USA.

M3533EN.2016-05 Rev. B
Printed 05/2016, XXXX