

A Focus on High Performance

HexaLux™ Exam Light




Brilliant. Powerful, Next-Generation LED Exam Light.

Experience the clear difference in coolness, comfort and natural color rendition inherent in pure, white LED light in your exam room. One light. Many possibilities.



Configured for Unlimited Clinical Applications

Its versatility in design includes five light intensity levels and configurations that mount to the ceiling, wall, rail or mobile unit. It's ideal for a multitude of rooms including exam room, birthing suite, or ICU.

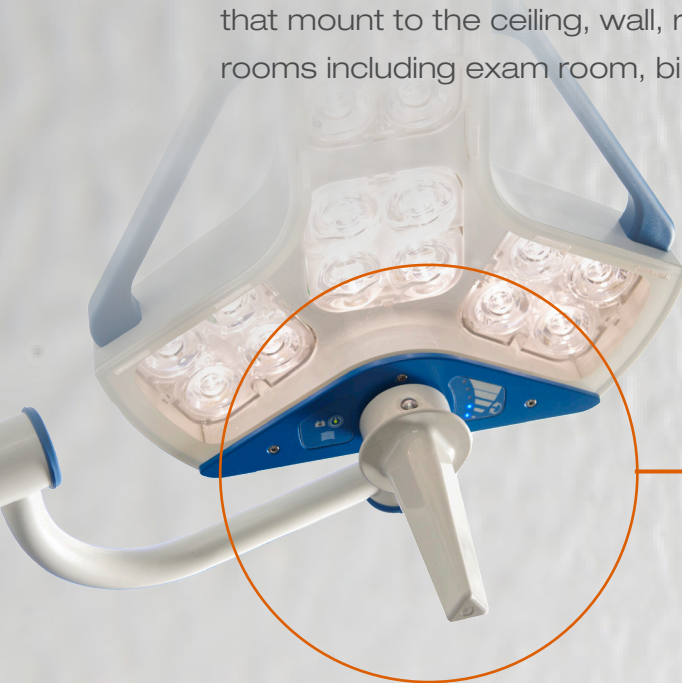

Sterilizable, Centrally Controlled Removable Handle

For simplified use by staff, all features are directly controlled from a central handle, which can be sterilized for use in minor surgical procedures.

- Turn light on or off
- Adjust light intensity
- Adjust light spot size



Further Customize Through Accessories

A convenient **shelf and basket** can be used for patient records, accessories, dressings, etc.

An optional **retractable power cord** eliminates cables on the floor.

Built to Last. And Last.

The trifecta has been accomplished; a **powerful**, high R9 value of 97, with a 65,000 Lux and long-lasting bulb life of up to 50,000 hours. With 16 LEDs its output provides an unmatched light source at a **competitive, affordable price**.

Easy Installation in Any Room Design

The modular design and partial assembly of the unit when it arrives, makes for quick installation and less room downtime.

Unlimited Clinical Applications

Flexibility in mounting and location within the room allow for Customer Adaptability



Ceiling Mount System



Wall Mount System

Technical Specifications

Features	HexaLux™ Exam Light
Light Intensity	Up to 65,000 Lux
Pattern Size	5.5" to 8.1" (14.0 to 20.5 cm)
Color Temperature	4,400 K (+/- 6%)
CRI (Color rendering index)	95 (+/- 6%)
R9	97 (+/- 6%)
Average LED Bulb Life	50,000 Hours

Mobile Stand

* Available with or without battery backup option



For more information, contact your STERIS representative or visit www.steris.com.



STERIS Corporation
5960 Heisley Road
Mentor, OH 44060-1834 · USA
800-548-4873
www.steris.com