

**APPLICATION**

STERIS offers a wide variety of material handling racks to provide flexibility and versatility for processing reusable materials in the following washer/disinfectors:

- Reliance® Vision® Single-Chamber Washer/Disinfectors
- Hamo® Vision® Single-Chamber Washer/Disinfectors<sup>1</sup>
- Reliance® Genfore™ Washer/Disinfectors
- Hamo® Genfore™ Washer/Disinfectors<sup>1</sup>
- Reliance® Vision® Multi-Chamber Washer/Disinfectors
- Hamo® Vision® Multi-Chamber Washer/Disinfectors<sup>1</sup>
- Reliance® Synergy® Washer/Disinfectors
- Hamo® Synergy® Washer/Disinfectors<sup>1</sup>
- Reliance® 444 Single-Chamber Washer/Disinfectors
- Reliance® 777 Automated Multi-Chamber Washer/Disinfectors

- AMSCO® 3052 Single-Chamber Washer/Disinfectors
- AMSCO® 5052 Single-Chamber Washer/Disinfectors

Diverse material processing requirements can be accommodated with these racks, and applications can be configured to the specific needs of the facility. Racks accommodate mixed loads, allowing a variety of instruments and utensils to be processed simultaneously.

**DESCRIPTION**

All accessories are constructed of stainless steel, are electropolished, and are not affected by detergents or by the heat and humidity produced during normal operation.

*NOTE: This Tech Data is divided into Current and Traditional Products.*

**Selections Checked Below Apply To This Equipment**

**ACCESSORIES FOR CURRENT PRODUCTS**

- (FD74-800) 2-Level Vision Manifold Rack  
Qty. \_\_\_\_
- (FD79-100) 2-Level Vision Manifold Rack (evenly split)  
Qty. \_\_\_\_
- (FD74-900) 3-Level Vision Manifold Rack  
Qty. \_\_\_\_
- (FD75-100) 4-Level Vision Manifold Rack  
Qty. \_\_\_\_
- (FD75-000) 4-Level Single-Chamber Manifold Rack  
Qty. \_\_\_\_
- (FDV72-000) 4-Level Suction Tips and Instrument Manifold Rack  
Qty. \_\_\_\_
- (FD75-200) 5-Level Vision Manifold Rack  
Qty. \_\_\_\_
- (FD77-900) 3-Level Manifold Rack with Removable MIS Upper Insert  
Qty. \_\_\_\_
- (FD02-6) 3-Level Manifold Rack (Increased Height) with Removable MIS Upper Insert  
Qty. \_\_\_\_
- (FD78-000) Vision Anesthesia/Respiratory Rack  
Qty. \_\_\_\_
- (FD75-500) Multi-Function Rack for Large Items  
Qty. \_\_\_\_
- (FD75-600) Multi-Function Rack for Small Items  
Qty. \_\_\_\_
- (FD77-100) Vision General Purpose Rack  
Qty. \_\_\_\_
- (CM53) General Purpose Rack  
Qty. \_\_\_\_
- (FD75-300) Flexible Hold-Down Screen  
Qty. \_\_\_\_
- (HA108340) Mesh Tray Lid  
Qty. \_\_\_\_

**ACCESSORIES FOR CURRENT PRODUCTS (Cont'd)**

- (247010029H) Instrument Mesh Tray  
Qty. \_\_\_\_
- (FD78-100) Instrument Stringers (20 Units for up to 10 Mesh Trays)  
Qty. \_\_\_\_
- (FD76-900) Flexible Hose for Cannulated Instruments  
Qty. \_\_\_\_
- (FD21-800) Transfer Cart  
Qty. \_\_\_\_
- (FD61-700) Universal Transfer Cart  
Qty. \_\_\_\_

**ACCESSORIES FOR TRADITIONAL PRODUCTS**

- (FD71-800) Rigid Hold-Down Screen  
Qty. \_\_\_\_
- (FD66-400) Rigid MIS Instrument Rack  
Qty. \_\_\_\_
- (CM54) Basin Rack  
Qty. \_\_\_\_
- (CM55) Tray Rack  
Qty. \_\_\_\_
- (CM56) Bedpan Rack with Spray Nozzles  
Qty. \_\_\_\_
- (CM57) Glassware Rack  
Qty. \_\_\_\_
- (FD70-400) Anesthesia/Respiratory Rack  
Qty. \_\_\_\_
- (FD24-300) Transfer Cart (Assembled)  
Qty. \_\_\_\_

**ACCESSORIES FOR PHASED-OUT WASHERS**

- (FDV70S-500) 2-Level Manifold Rack  
Qty. \_\_\_\_
- (FD72S-100) 2-Level Suction Tips and Instrument Rack, 777  
Qty. \_\_\_\_
- (FD57S-800) 2-Level Manifold Rack Assembly, 777  
Qty. \_\_\_\_
- (FD70S-500) 2-Level Manifold Rack Assembly, 444  
Qty. \_\_\_\_
- (FD70S-600) 3-Level Manifold Rack Assembly, 444  
Qty. \_\_\_\_
- (FDV24S-600) 4-Level Manifold Rack Assembly, 444  
Qty. \_\_\_\_
- (FDV24S-500) 5-Level Manifold Rack Assembly, 444  
Qty. \_\_\_\_
- (FD02S-000) Multiple-Purpose Basin Rack, 444/777  
Qty. \_\_\_\_
- (FD02S-100) Kidney Bowl Rack, 444  
Qty. \_\_\_\_
- (FD04S-000) Bedpan Rack, 444  
Qty. \_\_\_\_
- (FD60S-200) Transfer Cart, 777  
Qty. \_\_\_\_

*1. The Hamo line of Washers/Disinfectors, along with the registration to the Hamo trademark, is now owned by Amsonic AG, Switzerland.*

Item \_\_\_\_\_

Location(s) \_\_\_\_\_

\_\_\_\_\_

## CURRENT PRODUCTS MATRIX

**Table 1-1. Rack Compatibility Matrix for Current Products**

Equipment Numbers	Description	Washer/Disinfector			
		Reliance/Hamo Vision SC <sup>6</sup>	Reliance/Hamo Genfore/Synergy	Reliance/Hamo Vision MC <sup>6</sup>	AMSCO 3052 And 5052 SC <sup>6</sup>
<b>2-Level Manifold Racks</b>					
FD74-800	2-Level Vision Manifold Rack	Y	Y	Y	Y
FD79-100	2-Level Vision Manifold Rack (evenly split)	Y	Y	Y	Y
<b>3-Level Manifold Racks</b>					
FD74-900	3-Level Vision Manifold Rack	Y	Y	Y	Y
<b>4-Level Manifold Racks</b>					
FD75-100	4-Level Vision Manifold Rack	Y	Y	Y	Y
FD75-000	4-Level Single-Chamber Manifold Rack (Not compatible with Vision Multi-Chamber)	Y	Y	N	Y
FDV72-000	4-Level Suction Tips and Instrument Manifold Rack (See NOTE)	N	Y	Y	N
<b>5-Level Manifold Racks</b>					
FD75-200	5-Level Vision Manifold Rack	Y	Y	Y	Y
<b>Special Purpose Racks</b>					
FD77-900	3-Level Manifold Rack with Removable MIS Upper Insert	Y	N	N	N
FD02-6	3-Level Manifold Rack (Increased Height) with Removable MIS Upper Level	Y	N	N	N
FD78-000	Vision Anesthesia Respiratory Rack	Y	N	N	N
<b>General Purpose Racks</b>					
FD75-500	Multi-Function Rack for Large Items	Y <sup>1</sup>	Y <sup>1</sup>	Y <sup>1</sup>	Y
FD75-600	Multi-Function Rack for Small Items	Y <sup>2</sup>	Y <sup>2</sup>	Y <sup>2</sup>	Y
FD77-100	Vision General Purpose Rack	Y	Y	Y	Y
CM53	General Purpose Rack	Y <sup>1</sup>	Y	Y	Y
<b>Hold-Down Screens and Lids</b>					
FD75-300	Flexible Hold-Down Screen	Y <sup>5</sup>	Y	Y	Y
HA108340	Mesh Tray Lid	Y <sup>3</sup>	Y <sup>3</sup>	Y <sup>3</sup>	Y <sup>3</sup>
<b>Miscellaneous Items</b>					
247010029H	Instrument Mesh Tray	Y <sup>3</sup>	Y <sup>3</sup>	Y <sup>3</sup>	Y <sup>3</sup>
FD78-100	Instrument Stringers (20 Units for up to 10 Mesh Trays)	Y	Y	Y	Y
FD76-900	Flexible Hose for Cannulated Instruments	Y <sup>4</sup>	Y <sup>4</sup>	Y <sup>4</sup>	Y <sup>4</sup>
<b>Transfer Carts</b>					
FD21-800	Transfer Cart (Unassembled)	Y	Y	Y	Y
FD61-700	Universal Transfer Cart	Y	Y	Y	Y

**Notes:**

1. Must be used with FD74-800
2. Must be used with FD74-800 or FD74-900
3. Must be used with any multi-level rack
4. Must be used with FD74-900, FD75-100, FD75-000 or FD75-200
5. Must be used with FD75-600
6. SC= Single-Chamber; MC = Multi-Chamber

*NOTE: For FDV72-000 (4-Level Suction Tips and Instrument Manifold Rack) – Suction tips racks are for use only during instrument cycle processing. Sliding inlet in the washer/disinfector chamber automatically connects to the bottom of the rack and supplies wash solution/rinse water to the rotary spray arms and specialized spray nozzles. Spray arms are designed for thorough cleaning of the bottom and top load, while spray nozzles offer thorough cleaning of the cannulated suction tips. Luer lock basket is removable to accommodate different load capacities.*

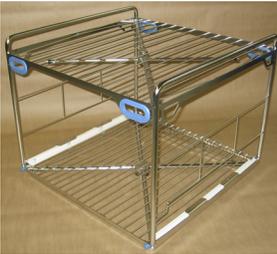
**ADVISORY NOTE:** STERIS does not recommend the use of this accessory for processing endoscopes, laparoscopes or other instruments used for minimally invasive surgery. Processing these types of instruments may result in equipment damage, as well as improperly cleaned/disinfected items.

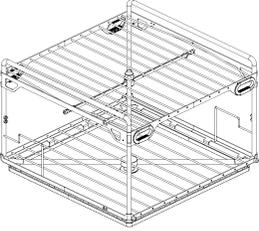
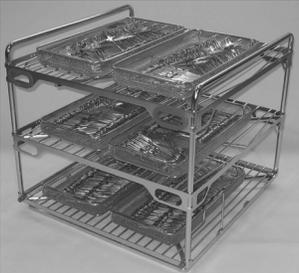
**Table 1-2. Rack Compatibility for Vision Multi-Level Racks**

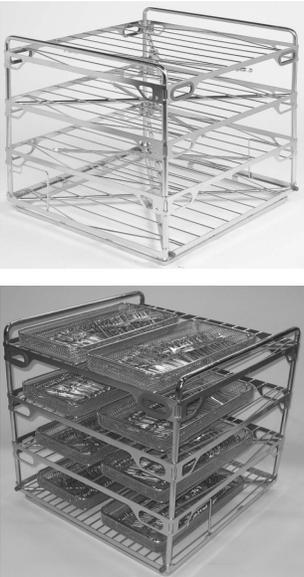
Multi-Level Compatibility Matrix for Vision Single-Chamber Washer		Vision Multi-Level Racks					
		2 Level		3 Level	4 Level	4 Level	5 Level
Equipment Number	Description	FD74-800	FD79-100	FD74-900	FD75-100	FD75-000	FD75-200
FD75-500	Multi-Function Rack for Large Items	Y	N	N	N	N	N
FD75-600	Multi-Function Rack for Small Items	Y	Y	Y	N	N	N
247010029H	Instrument Mesh Tray	Y	Y	Y	Y	Y	Y
FD76-900	Flexible Hose for Cannulated Instruments	N	Y	Y	Y	Y	Y

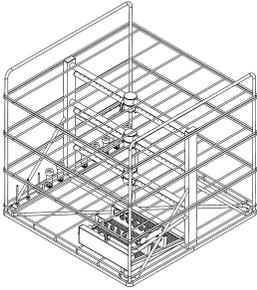
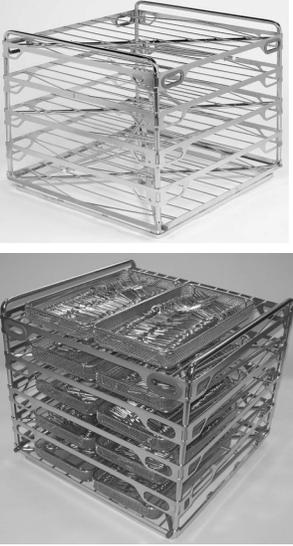
**MANIFOLD RACKS – CURRENT PRODUCTS**

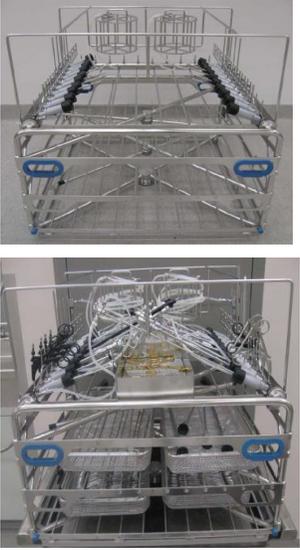
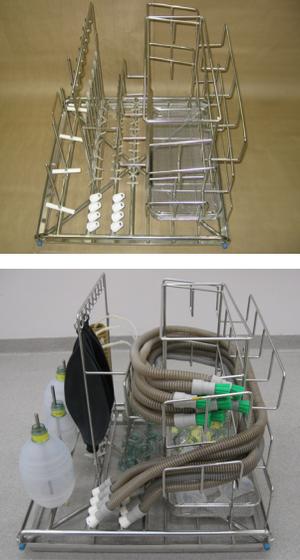
The manifold rack is used in conjunction with trays and racks not equipped with their own manifold. Sliding inlet in the washer/disinfector chamber automatically connects to the bottom of the manifold rack and supplies wash solution/rinse water to rotary spray arms or specialized spray nozzles. Spray arms are designed for thorough cleaning of the bottom and top of the load, while spray nozzles offer thorough cleaning of specific zones (e.g., inside of the anesthesia tubes).

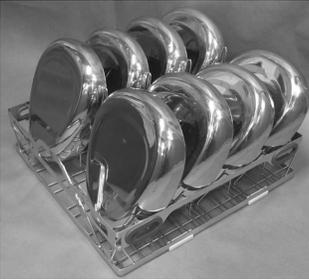
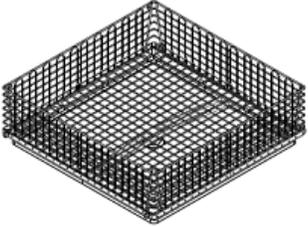
Accessory	Overall Dimensions W x L x H: Inches (mm) Weight: lb (kg)	Load Capacity Inches (mm)	Application					
<b>(FD74-800) 2-Level Vision Manifold Rack</b>  	25-1/2 x 25-1/2 x 22-1/4 (648 x 648 x 565)  25 (11)	<b>Top Level: (W x L x H)</b> 24-1/4 x 25-1/2 x see table below (616 x 648 x see table below)	Used to hold loaded instrument trays during processing.  Used to hold general purpose racks and special purpose racks during processing.  New Cool Touch™ Handles and Spray Arm design.  Accommodates: <ul style="list-style-type: none"> <li>• (FD75-500) Multi-Function Rack for Large Items</li> <li>• (FD75-600) Multi-Function Rack for Small Items</li> <li>• (FD75-300) Flexible Hold-Down Screen</li> <li>• (247010029H) Instrument Mesh Tray</li> <li>• (HA108340) Mesh Tray Lid</li> <li>• (FD71-800) Hold-Down Screen</li> </ul>					
		<table border="1"> <tr> <td></td> <td>Vision SC</td> <td>Vision MC</td> <td>AMSCO 3052 And 5052</td> </tr> <tr> <td>Top Level Clearance</td> <td>6-3/4 (171)</td> <td>4 (102)</td> <td>6-1/2 (165)</td> </tr> </table>			Vision SC	Vision MC	AMSCO 3052 And 5052	Top Level Clearance
	Vision SC	Vision MC	AMSCO 3052 And 5052					
Top Level Clearance	6-3/4 (171)	4 (102)	6-1/2 (165)					
		<b>Bottom Level: (W x L x H)</b> 24-1/4 x 25-1/2 x 15-3/4 (616 x 648 x 400)						

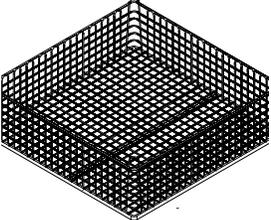
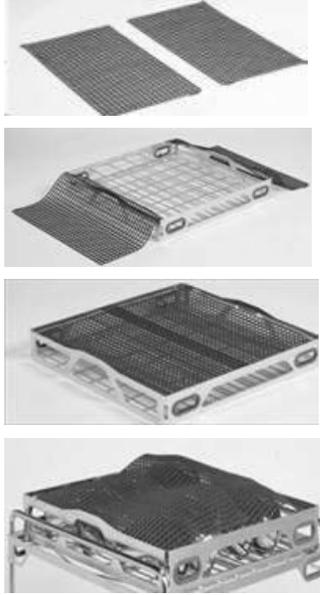
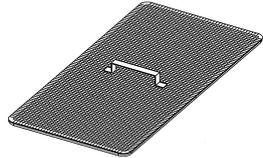
Accessory	Overall Dimensions W x L x H: Inches (mm) Weight: lb (kg)	Load Capacity Inches (mm)	Application								
<p><b>(FD79-100) 2-Level Vision Manifold Rack (evenly split)</b></p> 	<p>25-1/2 x 25-1/2 x 18 (648 x 648 x 457)</p> <p>25 (11)</p>	<p><b>Top Level:</b> (W x L x H) 24-1/4 x 25-1/2 x see table below (616 x 648 x see table below)</p> <table border="1" data-bbox="743 411 1201 575"> <tr> <td></td> <td>Vision SC</td> <td>Vision MC</td> <td>AMSCO 3052 And 5052</td> </tr> <tr> <td>Top Level Clearance</td> <td>11 (279)</td> <td>8-1/4 (209)</td> <td>10-3/4 (273)</td> </tr> </table> <p><b>Bottom Level:</b> (W x L x H) 24-1/4 x 25-1/2 x 11-1/2 (616 x 648 x 292)</p>		Vision SC	Vision MC	AMSCO 3052 And 5052	Top Level Clearance	11 (279)	8-1/4 (209)	10-3/4 (273)	<p>Used to hold loaded instrument trays during processing.</p> <p>Used to hold general purpose racks and special purpose racks during processing.</p> <p>New Cool Touch™ Handles and Spray Arm design.</p> <p>Accommodates:</p> <ul style="list-style-type: none"> <li>• (FD75-600) Multi-Function Rack for Small Items</li> <li>• (247010029H) Instrument Mesh Tray</li> <li>• (FD76-900) Flexible Hose for Cannulated Instruments</li> </ul>
	Vision SC	Vision MC	AMSCO 3052 And 5052								
Top Level Clearance	11 (279)	8-1/4 (209)	10-3/4 (273)								
<p><b>(FD74-900) 3-Level Vision Manifold Rack</b></p>  	<p>25-1/2 x 25-1/2 x 22-1/4 (648 x 648 x 565)</p> <p>31 (14)</p>	<p><b>Top Level:</b> (W x L x H) 24-1/4 x 25-1/2 x see table below (616 x 648 x see table below)</p> <table border="1" data-bbox="743 1050 1201 1213"> <tr> <td></td> <td>Vision SC</td> <td>Vision MC</td> <td>AMSCO 3052 And 5052</td> </tr> <tr> <td>Top Level Clearance</td> <td>7 (178)</td> <td>4-1/4 (108)</td> <td>6-3/4 (171)</td> </tr> </table> <p><b>Middle and Bottom Levels:</b> (W x L x H) 24-1/4 x 25-1/2 x 6-1/2 (616 x 648 x 165)</p>		Vision SC	Vision MC	AMSCO 3052 And 5052	Top Level Clearance	7 (178)	4-1/4 (108)	6-3/4 (171)	<p>Used to hold loaded instrument trays during processing.</p> <ul style="list-style-type: none"> <li>• Used to hold three (FD75-600) Multi-Function Racks for Small Items.</li> </ul> <p>Also holds:</p> <ul style="list-style-type: none"> <li>• (247010029H) Instrument Mesh Tray</li> <li>• (FD76-900) Flexible Hose for Cannulated Instruments</li> </ul> <p>New Cool Touch™ Handles and Spray Arm design.</p>
	Vision SC	Vision MC	AMSCO 3052 And 5052								
Top Level Clearance	7 (178)	4-1/4 (108)	6-3/4 (171)								

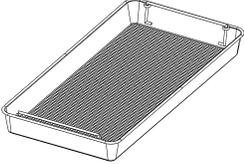
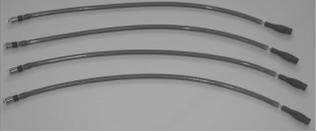
Accessory	Overall Dimensions W x L x H: Inches (mm) Weight: lb (kg)	Load Capacity Inches (mm)	Application								
<p><b>(FD75-100) 4-Level Vision Manifold Rack</b></p> 	<p>25-1/2 x 25-1/2 x 22-1/4 (648 x 648 x 565)</p> <p>39 (18)</p>	<p><b>Top Level:</b> (W x L x H) 24-1/4 x 25-1/2 x see table below (616 x 648 x see table below)</p> <table border="1" data-bbox="743 411 1201 575"> <tr> <td></td> <td>Vision SC</td> <td>Vision MC</td> <td>AMSCO 3052 And 5052</td> </tr> <tr> <td>Top Level Clearance</td> <td>6-3/4 (171)</td> <td>4 (102)</td> <td>6-1/2 (165)</td> </tr> </table> <p><b>Middle and Bottom Levels:</b> (W x L x H) 24-1/4 x 25-1/2 x 3-1/2 (616 x 648 x 89)</p>		Vision SC	Vision MC	AMSCO 3052 And 5052	Top Level Clearance	6-3/4 (171)	4 (102)	6-1/2 (165)	<p>Used to hold loaded instrument trays during processing.</p> <ul style="list-style-type: none"> <li>Used to hold four (FD75-600) Multi-Function Racks for Small Items.</li> </ul> <p>Also holds:</p> <ul style="list-style-type: none"> <li>(247010029H) Instrument Mesh Tray</li> <li>(FD76-900) Flexible Hose for Cannulated Instruments.</li> </ul> <p>New Cool Touch™ Handles and Spray Arm design.</p>
	Vision SC	Vision MC	AMSCO 3052 And 5052								
Top Level Clearance	6-3/4 (171)	4 (102)	6-1/2 (165)								
<p><b>(FD75-000) 4-Level Single-Chamber Manifold Rack (Not compatible with Vision Multi-Chamber Washer/Disinfector)</b></p> 	<p>25-1/2 x 25-1/2 x 23-3/4 (648 x 648 x 603)</p> <p>38 (17)</p>	<p><b>Top Level:</b> (W x L x H) 24-1/4 x 25-1/2 x see table below (616 x 648 x see table below)</p> <table border="1" data-bbox="743 1005 1198 1169"> <tr> <td></td> <td>Vision SC</td> <td>Vision MC</td> <td>AMSCO 3052 And 5052</td> </tr> <tr> <td>Top Level Clearance</td> <td>5-1/4 (133)</td> <td>N/A</td> <td>5 (127)</td> </tr> </table> <p><b>Middle and Bottom Levels:</b> (W x L x H) 24-1/4 x 25-1/2 x 4-1/4 (616 x 648 x 108)</p>		Vision SC	Vision MC	AMSCO 3052 And 5052	Top Level Clearance	5-1/4 (133)	N/A	5 (127)	<p>Used to hold loaded instrument trays during processing.</p> <ul style="list-style-type: none"> <li>Used to hold four (FD75-600) Multi-Function Racks for Small Items.</li> </ul> <p>Also holds:</p> <ul style="list-style-type: none"> <li>(FD76-900) Flexible Hose for Cannulated Instruments.</li> </ul> <p>New Cool Touch™ Handles and Spray Arm design.</p>
	Vision SC	Vision MC	AMSCO 3052 And 5052								
Top Level Clearance	5-1/4 (133)	N/A	5 (127)								

Accessory	Overall Dimensions W x L x H: Inches (mm) Weight: lb (kg)	Load Capacity Inches (mm)	Application								
<p><b>(FDV72-000) 4-Level Suction Tips and Instrument Manifold Rack</b></p> 	<p>24 x 24 x 23 (610 x 610 x 584) 41 (19)</p>	<p><b>Top Level:</b> 22-3/4 x 23-1/2 x 4-3/8 (578 x 597 x 111)</p> <p><b>Levels two and three:</b> With one item: 21-7/8 x 23-1/2 x 3-1/2 (556 x 597 x 89)</p> <p><b>Or, Two Instrument Trays:</b> 10-5/8 x 22-1/2 x 3-7/8 (270 x 572 x 98)</p> <p><b>Bottom Level:</b> With one item 21-1/4 x 23-1/2 x 3-1/2 (540 x 597 x 89)</p> <p><b>Or, Two Instrument Trays:</b> 10-5/8 x 21-1/2 x 3-7/8 (270 x 546 x 98)</p>	<p>Used to hold, during processing in the washer/disinfector:</p> <ul style="list-style-type: none"> <li>loaded instrument trays</li> <li>cannulated suction tips</li> <li>miscellaneous small items</li> </ul> <p>For use only during instrument cycle.</p> <p><b>Not compatible with Vision Single Chamber.</b></p>								
<p><b>(FD75-200) 5-Level Vision Manifold Rack</b></p> 	<p>25-1/2 x 25-1/2 x 22-1/4 (648 x 648 x 565)</p> <p>49 (22)</p>	<p><b>Top Level: (W x L x H)</b> 24-1/4 x 25-1/2 x see table below (616 x 648 x see table below)</p> <table border="1" data-bbox="743 842 1201 1010"> <tr> <td></td> <td>Vision SC</td> <td>Vision MC</td> <td>AMSCO 3052 And 5052</td> </tr> <tr> <td>Top Level Clearance</td> <td>6-1/4 (159)</td> <td>3-1/2 (89)</td> <td>6 (152)</td> </tr> </table> <p><b>Middle and Bottom Levels: (W x L x H)</b> 24-1/4 x 25-1/2 x 2-7/16 (616 x 648 x 62)</p>		Vision SC	Vision MC	AMSCO 3052 And 5052	Top Level Clearance	6-1/4 (159)	3-1/2 (89)	6 (152)	<p>Used to hold loaded instrument trays during processing.</p> <p>Used to hold:</p> <ul style="list-style-type: none"> <li>Five (FD75-600) Multi-function Racks for Small Items.</li> </ul> <p>Also holds:</p> <ul style="list-style-type: none"> <li>(247010029H) Instrument Mesh Tray and</li> <li>(FD76-900) Flexible Hose for Cannulated Instruments.</li> </ul> <p>New Cool Touch™ Handles and Spray Arm design.</p>
	Vision SC	Vision MC	AMSCO 3052 And 5052								
Top Level Clearance	6-1/4 (159)	3-1/2 (89)	6 (152)								

Accessory	Overall Dimensions W x L x H: Inches (mm) Weight: lb (kg)	Load Capacity Inches (mm)	Application
<p><b>(FD77-900)</b> 3-Level Manifold Rack with Removable MIS Upper Insert</p> <p><b>(FD02-6)</b> 3-Level Manifold Rack (Increased Height) with Removable MIS Upper Insert</p> 	<p>25-1/2 x 25-1/2 x 22 (647 x 647 x 558)</p> <p>49 (22)</p> <p>25-1/2 x 25-1/2 x 26 (647 x 647 x 660)</p> <p>49 (22)</p>	<p><b>Load Capacity Top Level:</b> Forty eight Cannulated instruments loaded on twenty instrument holders, eighteen large hollow jets and ten small hollow jets.</p> <p>Suction/irrigation ports can be connected with two flexible tubes for multi-port instruments and twelve flexible tubes for suction/irrigation connections.</p> <p>Two reels with ports are available for processing suction tubes, and one basket 7 x 23 x 2 (177 x 584 x 50) for reprocessing miscellaneous small items.</p> <p><b>Load Capacity Bottom Level:</b> Four Instrument Mesh Trays loaded with non-Cannulated instruments.</p>	<p>Used to process rigid minimally invasive surgical instruments (MIS).</p> <p>The top is used to hold Cannulated MIS instruments.</p> <p>The bottom is used to hold loaded instrument trays.</p>
<p><b>(FD78-000)</b> Vision Anesthesia/Respiratory Rack</p> 	<p>25-1/2 x 25-1/2 x 22 (647 x 647 x 558)</p> <p>29 (13)</p>	<p><b>Load Capacity:</b> Eight anesthesia tubes connected to eight Adapters for anesthesia tubing.</p> <p>Forty four anesthesia items and respiratory goods loaded on twelve large hollow jets, twelve medium hollow jets, twelve small hollow jets and eight hooks.</p> <p>One Mesh Tray loaded with small anesthesia items and respiratory goods.</p>	<p>Used to process different types of:</p> <ul style="list-style-type: none"> <li>• anesthesia tubing</li> <li>• bags</li> <li>• masks</li> <li>• respiratory goods</li> </ul>

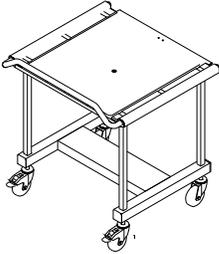
Accessory	Overall Dimensions W x L x H: Inches (mm) Weight: lb (kg)	Load Capacity Inches (mm)	Application								
<p><b>(FD75-500) Multi-Function Rack for Large Items</b></p>  	<p>23-1/2 x 23-1/2 x 12 (597 x 597 x 305)</p> <p>15 (7)</p>	<p><b>Top Level:</b> (W x L x H) 23 x 23 x see table below (584 x 584 x see table below)</p> <table border="1" data-bbox="743 411 1198 573"> <tr> <td></td> <td>Vision SC</td> <td>Vision MC</td> <td>AMSCO 3052 And 5052</td> </tr> <tr> <td>Top Level Clearance</td> <td>24-1/4 (616)</td> <td>21-1/4 (540)</td> <td>24 (610)</td> </tr> </table>		Vision SC	Vision MC	AMSCO 3052 And 5052	Top Level Clearance	24-1/4 (616)	21-1/4 (540)	24 (610)	<p>Used to hold:</p> <ul style="list-style-type: none"> <li>• basins</li> <li>• bowls</li> <li>• operating room shoes</li> <li>• bed pans, etc.</li> </ul> <p>Must be used with (FD74-800) 2- Level Vision Manifold Rack</p> <p>or</p> <p>(FD77-100) General Purpose Vision Rack.</p> <p>New Cool Touch™ Handles.</p>
	Vision SC	Vision MC	AMSCO 3052 And 5052								
Top Level Clearance	24-1/4 (616)	21-1/4 (540)	24 (610)								
<p><b>(FD75-600) Multi-Function Rack for Small Items</b></p>  	<p>23-1/2 x 23-1/2 x 4-1/2 (597 x 597 x 114)</p> <p>12 (6)</p>	<p><b>Top Level:</b> (W x L x H) 23 x 23 x see table below (584 x 584 x see table below)</p> <table border="1" data-bbox="743 1003 1198 1165"> <tr> <td></td> <td>Vision SC</td> <td>Vision MC</td> <td>AMSCO 3052 And 5052</td> </tr> <tr> <td>Top Level Clearance</td> <td>24-1/4 (616)</td> <td>21-1/4 (540)</td> <td>24 (610)</td> </tr> </table>		Vision SC	Vision MC	AMSCO 3052 And 5052	Top Level Clearance	24-1/4 (616)	21-1/4 (540)	24 (610)	<p>Used to hold:</p> <ul style="list-style-type: none"> <li>• small basins</li> <li>• bowls</li> <li>• light handles, etc.</li> </ul> <p>Must be used with:</p> <p>(FD74-800) or (FD74-900) Vision Multi-Level Rack</p> <p>or</p> <p>(FD77-100) General Purpose Vision Rack.</p> <p>New Cool Touch™ Handles.</p>
	Vision SC	Vision MC	AMSCO 3052 And 5052								
Top Level Clearance	24-1/4 (616)	21-1/4 (540)	24 (610)								
<p><b>(FD77-100) General Purpose Vision Rack</b></p> 	<p>25-1/2 x 25-1/2 x 8-1/2 (647 x 647 x 215)</p> <p>17 (7.7)</p>	<p><b>Top Level:</b> (W x L x H) 25 x 25 x see table below (635 x 635 x see table below)</p> <table border="1" data-bbox="743 1503 1198 1665"> <tr> <td></td> <td>Vision SC</td> <td>Vision MC</td> <td>AMSCO 3052 And 5052</td> </tr> <tr> <td>Top Level Clearance</td> <td>24 (609)</td> <td>21 (533)</td> <td>23-3/4 (603)</td> </tr> </table>		Vision SC	Vision MC	AMSCO 3052 And 5052	Top Level Clearance	24 (609)	21 (533)	23-3/4 (603)	<p>Used to hold two full size surgical instrument trays or miscellaneous hard goods.</p> <p>Used to hold items during processing in the washer/Disinfector.</p> <p>See note C on page 2.</p>
	Vision SC	Vision MC	AMSCO 3052 And 5052								
Top Level Clearance	24 (609)	21 (533)	23-3/4 (603)								

Accessory	Overall Dimensions W x L x H: Inches (mm) Weight: lb (kg)	Load Capacity Inches (mm)	Application								
<p><b>(CM53) General Purpose Rack</b></p> 	<p>24 x 24 x 7 (610 x 610 x 178)</p> <p>10 (5)</p>	<p><b>Top Level:</b> (W x L x H) 23-1/2 x 23-1/2 x see table below (597 x 597 x see table below)</p> <table border="1" data-bbox="743 411 1201 573"> <tr> <td></td> <td>Vision SC</td> <td>Vision MC</td> <td>AMSCO 3052 And 5052</td> </tr> <tr> <td>Top Level Clearance</td> <td>25 (635)</td> <td>22 (559)</td> <td>24-3/4 (629)</td> </tr> </table>		Vision SC	Vision MC	AMSCO 3052 And 5052	Top Level Clearance	25 (635)	22 (559)	24-3/4 (629)	<p>Used to hold items during processing in the washer/disinfector.</p> <p>Holds two full size surgical instrument trays or miscellaneous hard goods.</p> <p>Can be used alone in Genfore, Synergy or Vision Multi-Chamber washer/disinfectors.</p> <p>Must be used in (FD74-800) 2-Level Vision Manifold Rack with Vision single-chamber washer/disinfectors.</p>
	Vision SC	Vision MC	AMSCO 3052 And 5052								
Top Level Clearance	25 (635)	22 (559)	24-3/4 (629)								
<p><b>(FD75-300) Flexible Hold-Down Screen</b></p> 	<p>10 x 10 (254 x 254)</p> <p>3 (1)</p>	<p>N/A</p>	<p>Designed to be used with (FD75-600) Multi-Function Rack for Small Items to retain light items or serves as cushion for delicate items.</p> <p>Must be used with (FD74-800) 2-Level Vision Manifold Rack in Reliance/Hamo Vision Single-Chamber Washer/Disinfector.</p>								
<p><b>(HA108340) Mesh Tray Lid</b></p> 	<p>9-7/8 W x 17-1/2 L (251 x 445)</p> <p>1 (1/2)</p>	<p>N/A</p>	<p>Used with Instrument Mesh Tray to hold down lightweight plasticware.</p> <p>To be used with (247010029H) Instrument Mesh Tray.</p>								

Accessory	Overall Dimensions W x L x H: Inches (mm) Weight: lb (kg)	Load Capacity Inches (mm)	Application
<p><b>(227010029H) Instrument Mesh Tray</b></p> 	<p>9-7/8 x 19 x 2 (251 x 483 x 51)</p> <p>2 (1)</p>	<p>9-7/8 x 19 x clearance available into multi-level instrument manifold rack (251 x 483 x clearance available into multi-level instrument manifold rack)</p>	<p>Used to hold loaded instrument trays during processing in the washer/disinfector</p> <p>Must be used with a multi-level rack, e.g.:</p> <ul style="list-style-type: none"> <li>• FD74-800</li> <li>• FD74-900</li> <li>• FD75-100</li> <li>• FD75-200</li> </ul>
<p><b>(FD78-100) Instrument Stringers (20 Units for up to 10 Mesh Trays)</b></p>  	<p>9-7/8 x 19 x 2 (250 x 482 x 50)</p> <p>1 (1/2)</p>	<p>Fifty Hemostats</p>	<p>Used to hold instruments together during cycle and ensure better load coverage.</p> <p>Placed in instrument Mesh Tray.</p>
<p><b>(FD76-900) Flexible Hose for Cannulated Instruments (Four Hose Assemblies per Accessory)</b></p>  	<p>23 (584) long</p> <p>2 (1)</p>	<p>N/A</p>	<p>Used to process cannulated instruments.</p> <p>Can be installed only on 3, 4 and 5-level Vision manifold racks.</p> <p>Must be used with:</p> <ul style="list-style-type: none"> <li>• (FD74-900) 3-Level Vision Manifold Rack</li> <li>• (FD75-100) 4-Level Vision Manifold Rack</li> <li>• (FD75-200) 5-Level Vision Manifold Rack</li> </ul> <p>or</p> <ul style="list-style-type: none"> <li>• (FD75-000) 4-Level Vision Single-Chamber Manifold Rack</li> </ul>

## TRANSFER CARTS – CURRENT PRODUCTS

Stainless-steel carts are designed for transporting loaded accessories to and from the washer/disinfectors, as well as loading and unloading accessories. Carts can also be used to store trays and racks.

Accessory	Overall Dimensions W x L x H: Inches (mm) Weight: lb (kg)	Load Capacity Inches (mm)	Application
<p><b>(FD21-800) Transfer Cart</b> (Unassembled)</p> 	<p>29-7/8 x 28-3/4 x 40-1/4 (759 x 730 x 1022)</p> <p>44 (20)</p> <p>Loading Height: 31 (14)</p>	<p>One 24 x 24 (610 x 610) accessory</p>	<p>For transferring accessories into and out of the washer/disinfectors.</p> <p>If STERIS service is required, use SE10-405-20 for service assembly.</p>
<p><b>(FD61-700) Universal Transfer Cart</b></p> 	<p>29-3/4 x 30-3/4 x 29-1/4 to 37-1/2 (756 x 781 x 743 to 953)</p> <p>145 (66)</p> <p>Loading Height: 29-1/4 to 37-1/2 (743 to 953)</p>	<p>One 24 x 24 (610 x 610) accessory</p>	<p>Height-adjustable transfer cart for transferring accessories into and out of the washer/disinfectors.</p>

## TRADITIONAL PRODUCTS MATRIX

**Table 1-3. Rack Compatibility Matrix for Traditional Products**

Equipment Number	Description	Washer				
		Reliance/Hamo Vision Single-Chamber Washer/Disinfectors	Reliance/Hamo Genfore/Synergy Washer/Disinfectors	Reliance 444 Single-Chamber Washer/Disinfectors	Reliance/Hamo Vision Multi-Chamber Washer/Disinfectors	AMSCO 3052 And 5052 Single-Chamber Washer/Disinfectors
<b>2-Level Manifold Racks</b>						
FDV70-500 (NLA) (MTO# FDV70S-500)*	2-Level Manifold Rack	N	Y	Y	Y	N
FD72-100 (NLA) (MTO# FD72S-100)*	2-Level Suction Tips and Instrument Rack, 777 washer/disinfectors†	N	N	N	N	N
FD57-800 (NLA) (MTO# FD57S-800)*	2-Level Manifold Rack, 777 washer, disinfectors	N	N	N	N	N
FD70-500 (NLA) (MTO# FD70S-500)*	2-Level Manifold Rack	N	Y	Y	N	N
<b>3-Level Manifold Racks</b>						
FDV70-600 (NLA) (MTO# FD70S-600)*	3-Level Manifold Rack Assembly	N	Y	Y	Y	N
FD70-600 (NLA) (MTO# FD70S-600)*	3-Level Manifold Rack	N	Y	Y	N	N
<b>4-Level Manifold Racks</b>						
FDV24-600 (NLA) (MTO# FDV24S-600)*	4-Level Manifold Rack (without Central Post)	N	Y	Y	Y	N
FD24-600 (NLA)	4-Level Instrument Manifold Rack (without Central Post)	N	Y	Y	N	N
<b>5-Level Manifold Racks</b>						
FDV24-500 (NLA) (MTO# FDV24S-500)*	5-Level Manifold Instrument Rack	N	Y	Y	Y	N
FD24-500 (NLA)	5-Level Manifold Rack	N	Y	Y	N	N
<b>Hold-Down Screens and Lids</b>						
FD71-800	Rigid Hold-Down Screen	Y	Y	Y	Y	Y
<b>Special Purpose Racks</b>						
FD66-400	Rigid MIS Instrument Rack‡	N	Y	Y	N	Y
CM54	Basin Rack	YA	Y	Y	Y	Y
CM55	Tray Rack	YA	Y	Y	Y	Y
CM56	Bedpan Rack with Spray Nozzles	N	Y	Y	Y	Y
CM57	Glassware Rack	N	Y	Y	Y	N
FD70-400	Anesthesia/Respiratory Rack	N	Y	Y	Y	Y
FD02-000 (NLA) (MTO# FD02S-000)*	Multi-Purpose Basin Rack, 444/777 washer/disinfectors	YA	Y	Y	Y	N
FD02-100 (NLA) (MTO# FD02S-100)	Kidney Bowl Rack, 444/777 washer/disinfectors	YA	Y	Y	Y	N
FD04-000 (NLA) (MTO# FD04S-000)*	Bedpan Rack	Y	Y	Y	Y	N

**Table 1-3. Rack Compatibility Matrix for Traditional Products (Cont'd)**

Equipment Number	Description	Washer				
		Reliance/Hamo Vision Single-Chamber Washer/Disinfector	Reliance/Hamo Genfore/Synergy Washer/Disinfector	Reliance 444 Single-Chamber Washer/Disinfector	Reliance/Hamo Vision Multi-Chamber Washer/Disinfector	AMSCO 3052 And 5052 Single-Chamber Washer/Disinfector
<b>Transfer Carts</b>						
FD24-300 (NLA)	Transfer Cart Assembled	Y	Y	Y	Y	Y
FD60-200 (NLA) (MTO# FD60S-200)*	Transfer Cart TC for 777 washer/disinfector	N	N	N	N	N

<sup>A</sup> Must be used with FD74-800.

\* Made to Order (MTO)

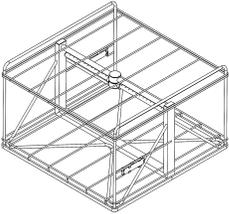
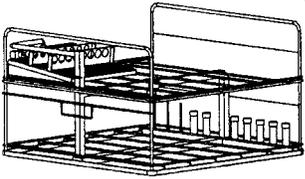
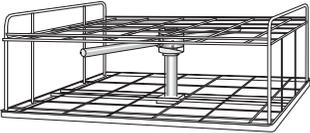
<sup>†</sup> Suction tips racks are for use only during instrument cycle processing. Sliding inlet in the washer/disinfector chamber automatically connects to the bottom of the rack and supplies wash solution/rinse water to the rotary spray arms and specialized spray nozzles. Spray arms are designed for thorough cleaning of the bottom and top load, while spray nozzles offer thorough cleaning of the cannulated suction tips. Luer lock basket is removable to accommodate different load capacities.

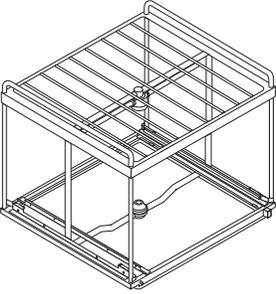
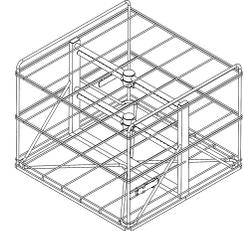
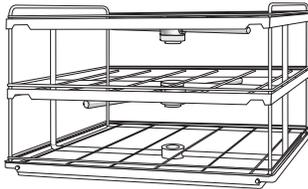
<sup>‡</sup> Rigid MIS Instrument Racks are for use only during Rigid MIS cycle processing. Sliding inlet in the washer/disinfector chamber automatically connects to bottom of rack and supplies wash solution/rinse water to rotary spray arms, specialized spraying jets and instrument holders.

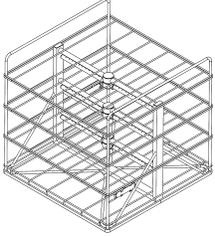
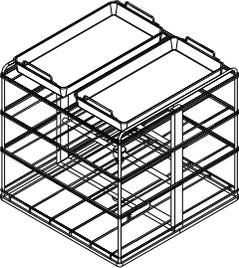
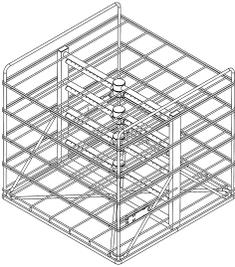
## MANIFOLD RACKS – TRADITIONAL PRODUCTS

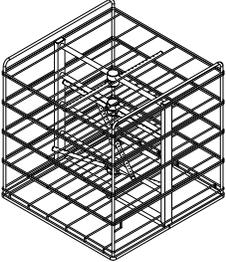
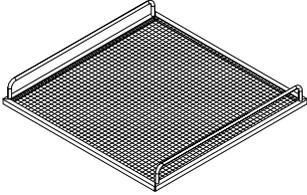
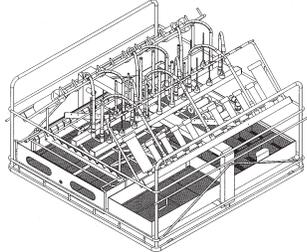
The manifold rack is used in conjunction with trays and racks not equipped with their own manifold. Sliding inlet in the washer/disinfector chamber automatically connects to the bottom of the manifold rack and supplies wash solution/rinse water to rotary spray arms or specialized spray nozzles. Spray arms are designed for thorough cleaning of the bottom and top of the load, while spray nozzles offer thorough cleaning of specific zones (e.g., inside of the anesthesia tubes).

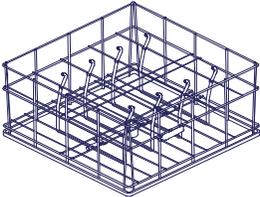
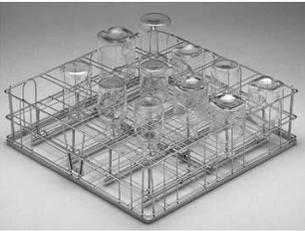
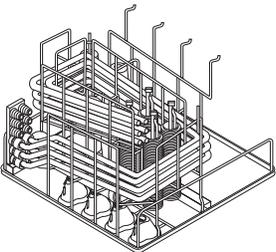
Made to Order (MTO) – approximately 45 days from order release.

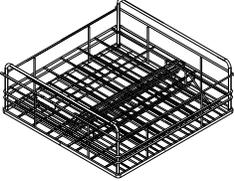
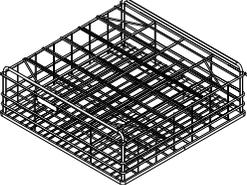
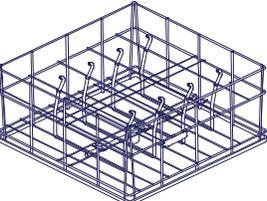
Accessory	Overall Dimensions W x L x H Inches (mm) Weight: lb (kg)	Load Capacity Inches (mm)	Application
<p><b>(FDV70-500) 2-Level Manifold Rack Assembly (NLA)</b> For 444 washer/disinfector, order MTO# FDV70S-500</p> <p>For all other washer/disinfectors, order FD74-800</p> 	<p>24 x 24 x 15 (610 x 610 x 381)</p> <p>22 (10)</p>	<p><b>Top Level:</b> (W x L x H) 22-3/4 x 23-1/2 x 10-1/4 (578 x 597 x 260)</p> <p><b>Lower Levels:</b> (W x L x H) <b>with One Item:</b> 21-7/8 x 23-1/2 x 9-1/8 (556 x 597 x 232)</p> <p><b>Or, Two Instrument Trays:</b> (W x L x H) 10-5/8 x 22-1/2 x 9-3/4 (270 x 572 x 248)</p>	<p>Used to hold, during processing in the washer/disinfector:</p> <ul style="list-style-type: none"> <li>• loaded instrument trays</li> <li>• general purpose baskets</li> <li>• bedpans</li> <li>• trays</li> <li>• basins or</li> <li>• glassware rack</li> </ul>
<p><b>(FD72-100) 2-Level Suction Tips and Instrument Rack (NLA)</b> For 777 washer/disinfector, order MTO# FD72S-100</p> 	<p>24 x 24 x 15 (610 x 610 x 381)</p> <p>22 (10)</p>	<p><b>Top:</b> Eight suction tips in luer lock basket and one full size tray 10-1/2 x 20-1/2 x 3 (267 x 521 x 76)</p> <p><b>Or:</b> one GPB</p> <p><b>Or:</b> two full size trays;</p> <p><b>Or:</b> one oversize tray 21 x 21 (533 x 533 max.)</p> <p><b>Bottom:</b> 10 suction tips on connectors and two full size trays;</p> <p><b>Or:</b> 10 suction tips on connectors and one oversize tray</p>	<p>Used to hold, during processing in the washer/disinfector:</p> <ul style="list-style-type: none"> <li>• loaded instrument trays</li> <li>• cannulated suction tips and</li> <li>• miscellaneous small items</li> </ul>
<p><b>(FD57-800) 2-Level Manifold Rack for 777 washer/disinfector (NLA)</b> For 777 washer/disinfector, order MTO# FD57S-800</p> 	<p>24 x 24 x 9 (610 x 610 x 229)</p> <p>20 (9)</p>	<p>Four full size trays 10-1/2 x 20-1/2 x 3 (267 x 521 x 76)</p> <p><b>Or:</b> Two full size trays and one oversize tray 22-1/2 x 22-1/2 x 4 (572 x 572 x 102), max.</p>	<p>Used to hold loaded instrument trays during processing in the washer/disinfector.</p>

Accessory	Overall Dimensions W x L x H Inches (mm) Weight: lb (kg)	Load Capacity Inches (mm)	Application
<p><b>(FD70-500) 2-Level Manifold Rack (NLA)</b> For 444 washer/disinfector, order MTO# FD70S-500</p> <p>For all other washer/disinfectors, order FD74-800</p> 	<p>26-1/2 x 25 x 22-1/8 (673 x 635 x 562)</p> <p>20 (9)</p>	<p><b>Top:</b> Two full size trays 10-1/2 x 20-1/2 x 3 (267 x 521 x 76) <b>Or:</b> One oversize tray 22-1/2 x 22-1/2 x 4 (572 x 572 x 102) max. <b>Bottom:</b> Any 24 x 24 (610 x 610) rack</p>	<p>Used to hold, during processing in the washer/disinfector:</p> <ul style="list-style-type: none"> <li>• loaded instrument trays</li> <li>• general purpose baskets</li> <li>• bedpans</li> <li>• trays and</li> <li>• basins</li> </ul> <p style="text-align: center;">or</p> <ul style="list-style-type: none"> <li>• glassware rack</li> </ul>
<p><b>(FDV70-600) 3-Level Manifold Rack Assembly (NLA)</b> For 444 washer, order MTO# FD70S-600</p> <p>For all other washer/disinfectors, order FD74-900</p> 	<p>24 x 24 x 19 (610 x 610 x 483)</p> <p>29 (13)</p>	<p><b>Top Level:</b> (W x L x H) 22-3/4 x 23-1/2 x 6-3/8 (578 x 597 x 162) <b>Lower Levels:</b> (W x L x H) <b>with One Item:</b> 27-7/8 x 23-1/2 x 5-1/2 (708 x 597 x 140) <b>Or: Two Instrument Trays:</b> (W x L x H) 10-5/8 x 22-1/2 x 5-7/8 (270 x 572 x 149)</p>	<p>Used to hold loaded instrument trays during processing in the washer/disinfector.</p>
<p><b>(FD70-600) 3-Level Manifold Rack (NLA)</b> For 444 washer, order MTO# FD70S-600</p> <p>For all other washer/disinfectors, order FD74-900</p> 	<p>26-1/2 x 25 x 22-1/8 (673 x 635 x 562)</p> <p>30 (14)</p>	<p><b>Top:</b> two full size trays (W x L x H) 10-1/2 x 20-1/2 x 6-3/4 (267 x 521 x 171) <b>Or:</b> one GPB <b>Or:</b> one oversize tray (W x L x H) 22-1/2 x 22-1/2 x 6-3/4 (572 x 572 x 171) max. <b>Middle:</b> two full size trays (W x L x H) 10-1/2 x 20-1/2 x 7-5/8 (267 x 521 x 194) <b>Or:</b> one GPB <b>Or:</b> one oversize tray (W x L x H) 22-1/2 x 22-1/2 x 7-5/8 (572 x 572 x 194) max. <b>Bottom:</b> two full size trays (W x L x H) 10-1/2 x 20-1/2 x 7-5/8 (267 x 521 x 194) <b>Or:</b> one GPB <b>Or:</b> one oversize tray (W x L x H) 22-1/2 x 22-1/2 x 7-5/8 (572 x 572 x 194) max.</p>	<p>Used to hold loaded instrument trays during processing in the washer/disinfector.</p>

Accessory	Overall Dimensions W x L x H Inches (mm) Weight: lb (kg)	Load Capacity Inches (mm)	Application
<p><b>(FDV24-600) 4-Level Manifold Rack (without central post) (NLA)</b> For 444 washer/disinfector, order MTO# FDV24S-600</p> <p>For Vision Single-Chamber washer/disinfector, order FD75-000</p> <p>For Genfore, Synergy, Vision Multi-Chamber washer/disinfectors, order FD75-100</p> 	<p>24 x 24 x 23 (610 x 610 x 584)</p> <p>34 (15)</p>	<p><b>Top Level:</b> (W x L x H) 22-3/4 x 23-1/2 x 4-3/8 (578 x 597 x 111) <b>Lower Levels:</b> (W x L x H) <b>with One Item:</b> 21-7/8 x 23-1/2 x 3-1/2 (555 x 597 x 89) <b>Or: Two Instrument Trays:</b> (W x L x H) 10-5/8 x 22-1/2 x 3-7/8 (270 x 572 x 98)</p>	<p>Used to hold loaded instrument trays during processing in the washer/disinfector.</p>
<p><b>(FD24-600) 4-Level Instrument Manifold Rack (without central post) (NLA)</b> For Vision Single-Chamber washer/disinfector, order FD75-000</p> <p>For Genfore, Synergy, Vision Multi-Chamber washer/disinfector, order FD75-100</p> 	<p>24 x 24 x 22-3/4 (610 x 610 x 578)</p> <p>35 (16)</p>	<p>Two full size instrument trays per level. Qty. 8</p>	<p>Used to hold loaded instrument trays during processing in the washer/disinfector.</p>
<p><b>(FDV24-500) 5-Level Manifold Instrument Rack Assembly (NLA)</b> For 444 washer/disinfector, order MTO# FDV24S-500</p> <p>For all other washer/disinfectors, order FD75-200</p> 	<p>24 x 24 x 23 (610 x 610 x 584)</p> <p>38 (17)</p>	<p><b>Top Level:</b> (W x L x H) 22-3/4 x 23-1/2 x 3-1/4 (578 x 597 x 83) <b>Lower Levels:</b> (W x L x H) <b>with One Item:</b> 21-7/8 x 23-1/2 x 2-3/8 (555 x 597 x 60) <b>Or: Two Instrument Trays</b> (W x L x H) 10-5/8 x 22-1/2 x 2-5/8 (270 x 572 x 67)</p>	<p>Used to hold loaded instrument trays during processing in the washer/disinfector.</p>

Accessory	Overall Dimensions W x L x H Inches (mm) Weight: lb (kg)	Load Capacity Inches (mm)	Application
<p><b>(FD24-500) 5-Level Manifold Rack (NLA)</b> For all washer/disinfectors, order FD75-200</p> 	<p>24 x 24 x 25 (610 x 610 x 635)</p> <p>38 (17)</p>	<p>Two full size instrument trays per level</p>	<p>Used to hold loaded instrument trays during processing.</p>
<p><b>(FD71-800) Rigid Hold Down Screen</b> Available for all washer/disinfectors</p> 	<p>24 x 24 x 1 (610 x 610 x 25)</p> <p>10 (5)</p>	<p>Cover for general purpose rack</p>	<p>Used to secure light objects, such as plastic items and similar materials, in a general purpose rack.</p>
<p><b>(FD66-400) Rigid MIS Instrument Rack</b> Available for Synergy, Genfore, 444, 3052 and 5052 washer/disinfectors Not compatible with Vision Single Chamber or Vision Multi-Chamber</p> 	<p>24 x 24 x 14-3/4 (610 x 610 x 375)</p> <p>38 (17)</p>	<p><b>Top:</b> 20 cannulated instrument holders, 12 flexible connectors for suction/irrigation connections and 20 spraying jets for cannulated instruments</p> <p><b>Bottom:</b> 12 cannulated instrument holders, four flexible connectors for suction/irrigation connections and one basket 7 x 23 x 2 (178 x 584 x 50) for miscellaneous small items. Bottom tray height clearance: 2" (50 mm)</p>	<p>Used to process rigid MIS instruments.</p>
<p><b>(CM54) Basin Rack</b> For all washer/disinfectors, order FD75-500</p> 	<p>24 x 24 x 11 (610 x 610 x 279)</p> <p>20 (9)</p>	<p>Each set supports and vertically separates two basins or bowls up to 12 (305) in diameter, <b>or</b> one basin or bowl up to 14-3/4 (375) in diameter</p>	<p>Used to hold items during processing in the washer/disinfectors.</p> <p>Must be used with FD74-800 in Vision Single-Chamber Washer/Disinfectors.</p>

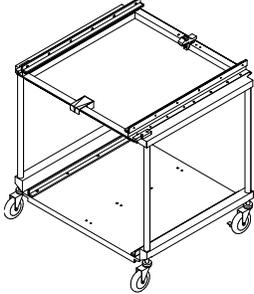
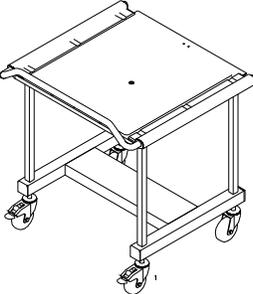
Accessory	Overall Dimensions W x L x H Inches (mm) Weight: lb (kg)	Load Capacity Inches (mm)	Application
<p><b>(CM55) Tray Rack</b> For all washer/disinfectors, order FD75-500 as replacement</p> 	<p>24 x 24 x 7 (610 x 610 x 178)  10 (5)</p>	<p>Accommodates nine 16 x 22 (406 x 559) <b>or</b> 18 11-1/2 x 15 (292 x 381) shallow trays</p>	<p>Used to hold items during processing in the washer/disinfecter.  Must be used with FD74-800 in Vision Single-Chamber Washer/Disinfecter.</p>
<p><b>(CM56) Bedpan Rack with Spray Nozzles</b> Available for Vision Multi-Chamber, Synergy, Genfore, 444, 3052 and 5052 washer/disinfectors. Not compatible with Vision Single Chamber.</p> 	<p>24 x 24 x 22-3/4 (610 x 610 x 578)  13 (5.9)</p>	<p>10-5/8 x 4-7/8 (270 x 123)</p>	<p>Used to hold up to eight bedpans upright during processing washer/disinfecter.</p>
<p><b>(CM57) Glassware Rack</b> For all washer/disinfectors, order FD75-600 as replacement</p> 	<p>24 x 24 x 7 (610 x 610 x 178)  12 (5)</p>	<p>16 sets of spacers to support an inverted beaker, bottle, jar, urinal or 75 to 2000 mL flask</p>	<p>Used to hold glassware during processing in the washer/disinfecter.</p>
<p><b>(FD70-400) Anesthesia/Respiratory Rack</b> Available for Vision Multi-Chamber, Synergy, Genfore, 444, 3052 and 5052 washer/disinfectors. Not compatible with Vision Single Chamber.</p> 	<p>24 x 24 x 20-3/4 (610 x 610 x 527)  30 (14)</p>	<p>12 tubes 3/4 to 1 dia. x 72 long (19 to 25 dia. x 1829 long) <b>Or</b> one tube 2-1/2 dia. x 19 long (64 dia. x 483 long), eight tubes 3/4 to 1 dia. x 72 long (19 to 25 dia. x 1829 long), four masks, four oxygen accumulators, four resuscitators, four ambulatory bags, four nebulizer bottles 7-1/4 (184) high, four nebulizer bottles 5 (127) high, one basket 10-1/2 x 22-1/2 x 3-1/2 (267 x 572 x 89)</p>	<p>Used to hold different types of tubing and medical devices during processing in the washer/disinfecter.</p>

Accessory	Overall Dimensions W x L x H Inches (mm) Weight: lb (kg)	Load Capacity Inches (mm)	Application
<p><b>(FD02-000) Multi-Purpose Basin Rack for 444/777 washer disinfectors (NLA)</b> For 444/777 washer/disinfectors, order MTO# FD02S-000</p> 	<p>24 x 24 x 9 (610 x 610 x 229)</p> <p>38 (17)</p>	<p>7-1/4 x 2 (184 x 51) qty. six 4-1/2 x 2 (114 x 51) qty. six 3-1/2 x 2 (89 x 51) qty. nine 7-1/2 x 5 (191 x 127) qty. six</p>	<p>Used to hold, during processing in the washer/disinfectors:</p> <ul style="list-style-type: none"> <li>• up to 27 basins of various shapes (round, square, rectangle) and sizes, as well as</li> <li>• kidney bowls</li> </ul> <p>Must be used with FD74-800 in Vision Single-Chamber Washer/Disinfectors.</p>
<p><b>(FD02-100) Kidney Bowl Rack for 444/777 Washer/Disinfectors (NLA)</b> For 444/777 washer/disinfectors, order MTO# FD02S-100</p> 	<p>24 x 24 x 8-3/8 (610 x 610 x 213)</p> <p>38 (17)</p>	<p>5-1/4 x 2-1/4 (133 x 57)</p>	<p>Used to hold up to 20 kidney bowls of various sizes, during processing in the washer/disinfectors.</p> <p>Must be used with FD74-800 in Vision Single-Chamber Washer/Disinfectors.</p>
<p><b>(FD04-000) Bedpan Rack without Spray Nozzles (NLA)</b> For 444 washer/disinfectors, order MTO#FD04S-000</p> <p>For all other washer/disinfectors, order FD75-500</p> 	<p>24 x 24 x 22-3/4 (610 x 610 x 578)</p> <p>12 (5.4)</p>	<p>10-5/8 x 4-7/8 (270 x 124)</p>	<p>Used to hold up to eight bedpans upright during processing in the washer/disinfectors.</p>

## **TRANSFER CARTS – TRADITIONAL PRODUCTS**

Stainless-steel carts are designed for transporting loaded accessories to and from the washer/disinfectors, as well as

loading and unloading accessories. Carts can also be used to store trays and racks.

Accessory	Overall Dimensions W x L x H Inches (mm) Weight lb (kg)	Load Capacity Inches (mm)	Application
<p><b>(FD24-300) Transfer Cart (Assembled)</b> Available for all washer/disinfectors</p> 	<p><b>MADE TO ORDER</b> 28-1/2 x 27 x 30-1/4 (724 x 686 x 768)</p> <p>50 (23)</p> <p>Loading Height: 31 (14)</p>	<p>Two racks <b>Or</b> one rack and one accessory header</p>	<p>For transferring accessories into and out of the washer/disinfectors. Includes one storage shelf.</p>
<p><b>(FD60-200) Transfer Cart, 777 Washer/Disinfectors (NLA)</b> For 777 washer/disinfectors, order MTO# FD60S-200</p> 	<p>29-7/8 x 28-3/4 x 40-1/4 (759 x 730 x 1022)</p> <p>44 (20)</p>	<p>One 24 x 24 (610 x 610) accessory</p>	<p>For transferring accessories into and out of washer/disinfectors.</p>

### **For Further Information, contact:**



STERIS Corporation  
5960 Heisley Road  
Mentor, OH 44060-1834 • USA  
440-354-2600 • 800-548-4873  
www.steris.com

*This document is intended for the exclusive use of STERIS Customers, including architects or designers. Reproduction in whole or in part by any party other than a Customer is prohibited.*